



NeoCrystal™ Sapphire Ceramic Bracket System

The "Crystal" Clear Advantage



THE “CRYSTAL” CLEAR ADVANTAGE

Manufactured from pure monocrystalline material, the **NeoCrystal** Sapphire Ceramic Bracket is the clear aesthetic choice for your practice. **NeoCrystal** Brackets are developed from a single crystal sapphire, one of the hardest and clearest materials found in nature. Once machined out of the crystal, the brackets are then annealed for strength and heat polished to perfect any microscopic flaws. The result is an unmatched aesthetic bracket that allows the natural tooth tone to show through, regardless of shade, for a virtually invisible appearance. Featuring **Henry Schein® Orthodontics™**’ proprietary **Crunch COAT™** Base and anatomically contoured base radii, the **NeoCrystal** Bracket provides the exceptional performance you demand.

UNMATCHED AESTHETICS

Pure monocrystalline material for a transparent, natural appearance

Stain resistant throughout treatment

DESIGNED FOR PATIENT COMFORT

Heat polished for a remarkably smooth surface

Low profile

Rounded corners

Single hook offers a less bulky feel than universal hooks




CRUNCH COAT BASE

Proprietary to *Henry Schein Orthodontics*, the Crunch Coat Base features:

Bond strength comparable to a mesh pad – the fine alumina particles bonded to the bracket base provide a large surface area for a strong adhesive bond

Consistent debonding – debonding occurs near the Crunch COAT adhesive interface providing minimal distress on the enamel surface

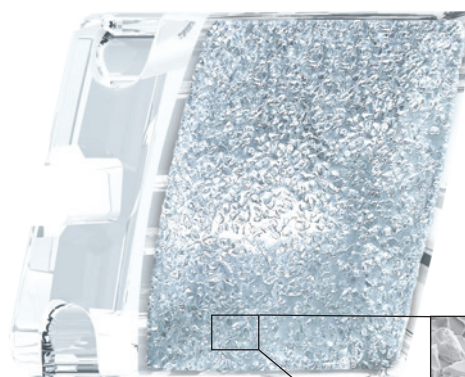

SEM Picture



ENHANCED PERFORMANCE

Heat treated for increased strength


Reduced friction offers sliding mechanics comparable to metal brackets

Non-allergenic for sensitive patients

Anatomically contoured base for quick and easy bracket placement


Removable color-code ID mark to ensure accurate placement

NEOCRystal BRACKETS ARE A PART OF THE NEOMAESTRO™ PREMIUM BRACKET SYSTEM-A FAMILY OF BRACKETS THAT ALL WORK IN UNISON.

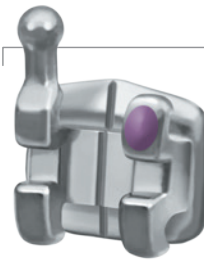

Maestro® Brackets for metal posterior bracketing—they offer reduced friction, improved control, and patient comfort.

NeoLucent® Plus™ Brackets for aesthetic anterior bracketing—they offer an invisible appearance, predictable debonding, and a low profile.

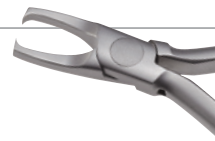

Triumph™ Debonding Pliers for metal posterior bracketing—they offer reduced friction, improved control, and patient comfort.

Triumph Angulated Debonding Pliers. The angulated design quickly and easily removes anterior and posterior direct bond brackets.



The mission of *Henry Schein Orthodontics* is to provide state-of-the-art orthodontic products and innovative clinical solutions that enable our customers to offer exceptional patient care, while expanding the scope and profitability of their practices.



For more information on our products and educational offerings, please contact us:
HenryScheinOrtho.com
Customer Service: 888.851.0533
In the U.S.: 800.547.2000
Outside the U.S.: +(1) 760 448 8600
Canada: CERUM 800.661.9567

To email an order:
USASales@HenryScheinOrtho.com
To email an order outside the U.S.:
IntlOrders@HenryScheinOrtho.com



 **HENRY SCHEIN®**
ORTHODONTICS

The New Movement In Orthodontics™

Ortho Organizers, Inc.
1822 Aston Ave.
Carlsbad, CA 92008-7306 USA