



## Maestro® Low-Profile Bracket System

Orchestrating controlled treatment  
for finely tuned finishes.



# Discover the Precision, Comfort, and Control of MIM-Manufactured *Maestro* Brackets

The *Maestro* Low-Profile Bracket System harmonizes the enduring precepts of classic orthodontic design with state-of-the-art engineering. The result is a complete orthodontic appliance system with corresponding buccal tubes—contemporary in form and function—that achieves the optimum clinical outcome in metal aesthetics, as well as the comfort that patients seek without sacrificing the precision and control that practitioners need. A finely tuned appliance design for predictable, efficient treatment, and discerning finishes.

## THE *MAESTRO* BRACKET OFFERS MULTIPLE VISUAL CUES TO ENHANCE PRECISE BRACKET PLACEMENT

The geometries of the *Maestro* Bracket are specifically designed to provide signposts which assist you in positioning brackets in their ideal location.

**Rhomboid-shaped** bracket and pad

**Lowest possible profile** reduces occlusal interferences

**Long-axis scribe line** not only runs the entire length of the pad, but also the bracket for an accurate placement cue



**Constant occlusal-gingival heights** across the system foster placement repeatability

**Vertical aspects** follow long axis of the tooth

**Horizontal aspects** run parallel to occlusal plane



**Specific bracket geometries**, such as in this lower lateral bracket, further guide the eye for placement accuracy

**Contoured archwire slot edges mesial-distal** reduce friction

**In/outs built-in, not machined** for precision and consistency

**Small footprint** is less bulky than conventional twin brackets, minimizing the look and feel



## IN-DEPTH RESEARCH ON TOOTH MORPHOLOGY SUPPORTS SUPERIOR ENGINEERING FOR A FINELY TUNED DESIGN AND FINISHES

**Smooth rounded surfaces** ensure patient comfort



**Anatomically contoured base and pad** nestle each tooth for positioning that feels right

Small profile is balanced with **deep, tie-wing undercuts** that ease secure ligation



Optional **gingivally offset brackets, on large bonding bases** of bicuspid, minimize bond failures

**80-gauge, micro-etched bondable mesh pad** offers proven bond adhesion



**Torque-in-base** provides level slot lineup, and helps avoid occlusal interferences


**Optimized mesial-distal width** offers superb rotational control



The *Maestro* Bracket is available in the Roth prescription, and the McLaughlin, Bennett, and Trevisi prescription, in .018" and .022" lumen sizes. Ball hooks are optional on laterals, cuspids, and bicuspid brackets; and each bracket is color-coded with an ID mark for easy identification.

## MAESTRO SYSTEM: A COMPLETE 7X7 SYSTEM

*Maestro* Buccal Tubes complete the appliance system. The 80-gauge mesh pad, compound contoured base, and defined buccal indent, provides you with the precise control needed throughout all phases of treatment. The large funneled entrance, side-grip areas, and color-coding make for easier identification, placement, positioning, plus wire insertion. The patient will have one of the most comfortable experiences with its smooth shape, low profile, rounded edges, and lowered hook. The *Maestro* Buccal Tube orchestrates all components of the system seamlessly, providing reduced chair time and increased satisfaction for both practitioner and patient.

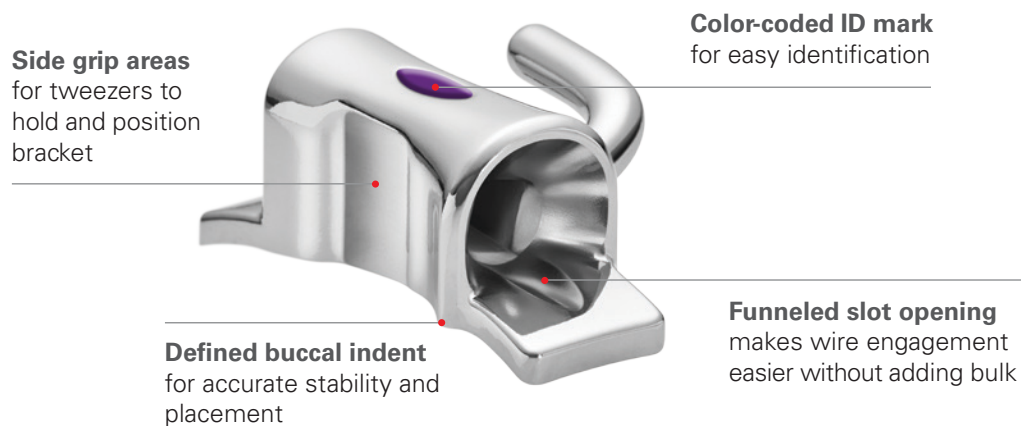

***Maestro Buccal Tubes*** are engineered to provide  
*a complete and seamless appliance system.*



For more information on our products and educational offerings, please contact us:  
HenryScheinOrtho.com  
Customer Service: 888.851.0533  
In the U.S.: 800.547.2000  
Outside the U.S.: +(1) 760 448 8600  
Canada: CERUM 800.661.9567

To email an order:  
USASales@HenryScheinOrtho.com  
To email an order outside the U.S.:  
IntlOrders@HenryScheinOrtho.com



**HENRY SCHEIN®**  
ORTHODONTICS

*The New Movement In Orthodontics™*

Ortho Organizers, Inc.  
1822 Aston Ave.  
Carlsbad, CA 92008-7306 USA